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declares that the object of the declaration is the product:

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conforms to essential requirement according to ECC Directive EMC 2014/30/UE, LVD 2014/35/UE, ROHS 2011/65/UE and DECRETO 15 October 2015, n. 206 applicable to this product. In addition, the object of the declaration described above, is in conformity with the relevant Union harmonization legislation:

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This declaration of conformity is issued under the sole responsibility of the manufacturer.

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Additional information:

Torino, 29/04/2016

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Signed for and on behalf of:

Walter Bottega, Importatore

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ENGLISH

Dear customer,

thank you for choosing Pressoestetica® Mesis®. You will appreciate the high quality materials and features that set this model apart.

We hope that you find this machine comfortable to use and remind you to read the manual, the contraindications and warnings prior to use.

The product must be used in accordance with the instructions contained in this manual.

The manufacturer and distributor do not accept any liability for any use that differs from the one indicated in this manual.

PROFESSIONAL USE: professionals who purchase and use equipment for professional purposes are directly liable both legally and criminally. In particular, please note that professional use of equipment by aestheticians is regulated by Ministerial Decree-Law 1/90.

Laws and regulations, to which reference is made in this manual, are European and Italian.

VIMA di W. B.
Via Pigafetta 65, 10129 Torino
fax 011.07.04.904
email: info@mesis.eu
www.mesis.eu

Reproduction, even in part, of any photo or text contained within this manual is prohibited.

© Copyright VIMA DI W.B.

Index

Index

Warnings	26
What is Pressoestetica®	34
Pressoestetica® Mesis® JoySense®2.0 - mod: JS02.....	36
Pressoestetica® Mesis® JoySense® 3.0 - mod: JS03.....	37
Pressoestetica® Mesis® JoySense® instructions for Use	38
Maintenance and Troubleshooting	41
Warranty	42
EU declaration of conformity for aesthetic use	43

 **Warnings**

Before using the device, please read the instructions for use and, in particular, the safety notes and follow these at all times.

Keep this manual with the device for future consultation for the entire product life. In the event that the product is sold to a third party, all documentation should also be handed over.

After removing the device from its packaging, check its integrity with the drawing and any presence of damage as a result of transportation. If in doubt, do not use the device and contact an authorised customer service centre.

Do not let children play with the packaging! Keep the plastic bag out of the reach of children; risk of suffocation!

Before connecting the device, please check that the mains voltage indicated on the identification label corresponds to those of the power supply available.

The identification label is located on the device.

This device must only be used for the purpose for which it was designed, i.e. as a machine for presotherapy. Any other use shall be deemed as non-compliant and therefore dangerous.

In order to ensure the best possible protection, it is advised to install a device with a differential current into the electric circuit, which has a nominal operating differential current not exceeding 30mA (ideally 10mA). Ask your installer for advice.

WARNING: do not use this device in the vicinity of the bathtub, sinks or other recipients containing water! Never immerse the device in water!

This product is not intended for use by persons (including children) with reduced physical, mental or sensory skills, or lack of experience and knowledge, unless there is the supervision or training on the use of the device supplied by the person responsible for their safety.

DO NOT use the device with wet hands or damp or bare feet.

DO NOT pull the power cord or the device itself when removing the plug from the socket.

DO NOT expose the device to humidity or to the influence of weather conditions (rain, sun).

DO NOT block the air holes on the device and do not place the device on a soft surface such as a bed or a sofa where the air holes could be blocked.

Do not cover the air holes.

Always pull the plug from the socket before cleaning or performing maintenance and when not in use.

In case of failure or malfunction of the device, switch it off and do not tamper with it. For any repairs, contact an authorised technical support centre only.

If the power cord is damaged, it must be replaced by an authorised technical support centre in order

to prevent any hazard.

To clean the device, use a dry cloth and check as present in this manual.

Please note that the device must only be used indoors and in accordance with the operating instructions. If the power cord is faulty, it must be replaced by a VIMA di B.W. customer service centre.

The electrical system to which the device is connected must conform with Ministerial Decree 37/08 and all safety measures in the system must be provided for such as the thermomagnetic circuit breakers on the main power lines.

Store and handle the product carefully. Use of the product is prohibited if it is damaged.

Do not place near sources of heat.

Do not use strong or abrasive detergents to clean the surface; use a dry cloth to clean the device when necessary.

If it is difficult to remove any dirt, use a mild detergent and wipe gently.

Do not use on patients with a pacemaker or internal defibrillator.

The parts that come into contact with the skin must be cleaned between treatments using the usual cleaning methods with mild disinfectant on a soft cloth.

Disconnect from the power supply before cleaning.

All parts that come into contact with the patient during treatment should be suitably sterilised and/or disinfected.

Never use the product on swollen or inflamed areas or damaged skin.

Use of the product is limited to residential/small office type environments.

Photos and features of the product are accurate at the time of printing, minor differences may be likely during production. Any differences do not affect the safety or performance of the product.

Warning:

To avoid the risk of fire or electric shock of any kind, do not expose the device to the rain or liquids of any kind.

The manufacturer shall provide information on the proper use of items. You should therefore always follow these instructions to ensure the product is properly and correctly installed and operated; as such they should be kept for future reference.

- Turn the switch to 0 (zero) or OFF (to turn the device off, press the **POWER** button) after use.
- Clean items with a dry cloth and do not use solvents or abrasive products. Avoid spilling liquids on the electrical parts.
- All of the device's technical specifications must be borne in mind.
- Check the symbols located on the identification label of each item and the icons visible in the (page [30-31-32-33](#) of this user manual) safety instruction section.
- Disposal of old products: the product has been designed and manufactured with high quality materials and components that can be recycled and reused. A crossed-through wheely bin symbol on your product indicates that the product meets the requirements of Community Directive WEEE 2012/19/EU.

Find out about the collection methods for electrical and electronic products in force in the area in which you wish to dispose of the product.

Please follow local regulations for the disposal of waste and do not dispose of old products with normal household waste. The proper disposal of products helps to prevent environmental pollution and possible health hazards.

This product is for indoor use only.

It is not suitable for installation in bathrooms.

The item must be connected directly to the mains through the appropriate power supply unit only.

Connection Y: in case of damage, the cord must only be replaced by the manufacturer, the distributor or by a technician to avoid any risks.

During installation ensure you do not make holes near the existing electrical wiring in the wall or ceiling.

Connecting cables must never be twisted and squashed into the space between the device and the supporting surface.

This product is not suitable for children.

Do not treat:

- carriers of active implants such as pacemakers, insulin pumps or other devices;
- carriers of metallic prostheses;
- subjects with inflammatory processes in place, with skin lesions, with neoplasms;
- pregnant women;
- do not treat subjects with sensitive skin, who may be restricted to powers of very low disbursement;
- Apply only on undamaged skin.

Operating mode:

The application time varies according to the treatment to be performed and, as a rule, can vary between 15 and 30 minutes.

In addition to the indications above, carefully follow the instructions, warnings and precautions for use described in the user manual.

The pressoestetica® JoySense® machine is designed, manufactured and tested in accordance with the strictest European safety standard regulations.

Any damage caused by use of the device in exceptional circumstances (coastal areas, industrial areas, frequent contact with organic fertiliser etc.) is not covered by the warranty conditions.

The warranty period starts from the date of purchase and is valid only upon presentation of an invoice, tax receipt or distance-purchase sales document (online shops, mail order catalogues). The warranty is void if the item was not installed according to the required instructions, or if repaired or tampered with by unauthorised personnel.

The manufacturer is also not liable for damage caused by improper use or incorrect application of the instructions.

Alterations in the design or technical specifications are reserved.

Environmental warnings

At the end of the product's useful life, dispose of it according to local regulations in force.

- if you use the product in a bathroom (STRONGLY DISCOURAGED), you must unplug it after use since proximity to water can be dangerous even when the product is switched off;

- In order to ensure the best possible protection, it is advised to install a device with a differential current into the electric circuit supplying the bathroom, which has a nominal operating differential current not exceeding 30 mA (ideally 10mA). Ask your installer for advice.

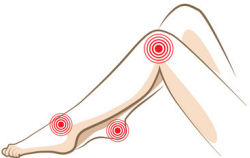
Warnings

- Before use, carefully read the manual in its entirety;
- Keep this manual close to the pressoestetica® device, so you can consult it at any time;
- This manual explains in detail how to properly use the device in order to prevent danger and damage to persons or property.

Meaning of the symbols in this manual

 Warning	 Precautions
<p>Failure to follow the instructions can result in serious injury or loss of life.</p>	<p>Failure to follow the instructions can result in the user causing injury to him/herself or damage to the device.</p>
Indications	
 <p>Indicates a restriction</p>	 <p>Is an indication to be followed properly</p>
 <p>Prohibits disassembly</p>	 <p>Indicates to unplug the power cord from the wall outlet</p>

 Warning	
<p>Do not use the product in the following cases:</p>  <p>People suffering from both internal and external inflammation or wounds at risk of bleeding must not use the product including those suffering from purulent inflammation. Do not use in case of general malaise, flu or illness.</p>	<p>Contraindications</p> <ul style="list-style-type: none"> • Pulmonary oedema • Thrombophlebitis • Heart failure • Infections • Deep vein thrombosis • Pulmonary embolism episodes • Wounds, lesions or tumours in the area of application or in proximity to areas where increased venous and lymphatic return could be harmful • Carriers of active implants such as pacemakers, insulin pumps or other devices; • Carriers of metallic prostheses; • Subjects with inflammatory processes in place, with skin lesions, with neoplasms; • Pregnant women.

	<p>Consult your doctor to check whether you have any disease which is incompatible with the aesthetic pressotherapy treatment. Any incompatible diseases are listed in the „<u>Contraindications section</u>“.</p>	
---	---	---

Warnings

Warning

	<p>If you feel unwell or experience pain during use, do not use the device. Discontinue use immediately.</p>		<p>Do not use the device for more than 15 minutes during the first treatment.</p> <p>Gradually increase the time in subsequent sessions depending on how well you tolerate the treatment. Begin the treatment starting with minimal pressure.</p>
	<p>Use for children is prohibited. Keep out of the reach of children. Risk of electric shock or other damage.</p>		<p>Do not use water or other substances that could come into contact with internal parts of the device.</p>
	<p>In case of failure or malfunction, do not open or attempt to repair the device. Risk of electric shocks and injuries to persons or property.</p>		<p>Your hands must be dry when connecting or disconnecting the power cord from the power supply (socket). Risk of electric shocks and injuries to persons or property. Risk of fire.</p>
	<p>Use only single-phase 220V current. After using the device, ensure you disconnect the power cord from the power supply properly. Risk of electric shocks and injuries to persons or property. Risk of fire.</p>		<p>Leave the power cord free, do not place it on top of, squash it against or attach it to any object. Risk of shock or electrical leakage. Risk of fire.</p>

Warnings

Warning

	<p>Use a plug socket for the exclusive use of the device. Too many plugs could cause the start of a fire, electrical leakage and electric shocks.</p>		<p>Do not hit the device or let it fall. This could cause serious damage to the device.</p>
	<p>Do not insert any object or product into the device. Risk of damage to the device and to people.</p>		<p>Do not place the device and its accessories next to sources of heat. It could cause serious damage to the device and/or accessories</p>
	<p>When wearing the leggings, abdominal or armband, do not move or walk.</p>		<p>Do not place the device and its accessories in direct sunlight or near the stove (gas, electric or any other fuel). This may cause deterioration, deformation and discoloration of the product.</p>
	<p>Do not use the device for a purpose other than massage. Risk of damage to people and the device.</p>		<p>Do not bend, knot, cut or break the power cord.</p>

Warnings

Warning

	<p>Do not place any object on top of the device. Do not place weights on the machine or on any of its connections, nor on electric cables and connecting air tubes. It may cause changes in the pressure.</p>		<p>If the machine has not been used for a while or is used in a cold place, it may not work properly, therefore, once switched on, leave to work at a low pressure for at least 3 inflation cycles. Do not leave the machine in a damp place or in the rain.</p>
	<p>Do not use water, oil, gasoline or chemicals to clean the accessories (leggings, armbands, etc.). Use a damp cloth. It could damage them, shortening their lifespan.</p>		<p>Do not use the accessories directly on the skin. Wear suitable undergarments underneath the applicators.</p>
	<p>Do not use water, oil, gasoline or chemicals to clean the accessories. Use a damp cloth.</p>		<p>Do not bend, knot, or cut the air tubes.</p>

Product subject to current provisions on waste sorting (WEEE Directive). Dispose of in the designated areas.

What is Pressoestetica®

The body requires daily attention, Mesis® provides you with effective, stylish and beneficial equipment to pamper you in the right way by strictly following the indications described in the user manual.

Pressoestetica® Mesis® JoySense® derives straight from the aesthetic pressotherapy systems.

In fact it corresponds to a treatment (equivalent to a manual massage) of aesthetic lymphatic drainage that acts locally, to counteract cellulite and water retention.

Aesthetic applications

- Adjuvant slimming
- Cellulite
- Stretch marks
- Aesthetic lymphatic drainage leg massage
- Aesthetic lymphatic drainage arm massage
- Water retention
- Heavy legs
- Swollen feet and legs
- Relaxation

Do not treat:

- carriers of active implants such as: pacemakers, insulin pumps or other devices;
- carriers of metallic prostheses;
- subjects with inflammatory processes in place, with skin lesions, with neoplasms;
- pregnant women;
- do not treat subjects with sensitive skin, who may be restricted to powers of very low disbursement;
- apply only on undamaged skin.

Main contraindications

Pulmonary oedema. Thrombophlebitis. Heart failure. Infections. Deep vein thrombosis. Pulmonary embolism episodes. Wounds, lesions or tumours in the area of application or in proximity to areas where increased venous and lymphatic return could be harmful. Carriers of pacemakers. Lymphangitis. Plexopathy and neuropathy.

Consult your doctor to check whether you have any disease which is incompatible with the aesthetic pressotherapy treatment. Any incompatible diseases are listed in the „Contraindications section“ of page 30 in sections „Do not treat“ and „Main contraindications“ of this page.

Index

Features and terminology for pressoestetica® Mesis® JoySense®:

- The pressoestetica® Mesis® JoySense® are **sequential** and intermittent;
- **Program** - chamber inflation mode;
- **Chambers or Sectors** - elements that comprise the applicator and that fill with air, in the pressoestetica® Mesis® JoySense® there are 5 chambers for every applicator;
- **Cycle** – one cycle is equivalent to the inflation of all chambers, from the first to the last regardless of the program used. After the last chamber has inflated, there is a pause and then a new cycle starts;
- **Applicator** – pneumo relaxation element composed of several chambers that inflate sequentially regulated by the chosen program (leggings, armband and abdominal/buttocks band);
- **Time** - duration of the session (adjustable by 15 – 25 - 30 minutes);
- **Pressure** - the quantity of air blown in the chamber to perform lymphatic drainage on the person, the pressure is adjustable and continues on each chamber.
- **Connector** - connecting tubes between the machine and the applicators. Thanks to “one touch connect”, inserting the connecting tubes for the air to the machine is easy and quick.

The operating cycle of a pressure exerted through use of the pressotherapy **follows the physiological movement of the blood and lymph** in a distal-proximal way. From the extremities to the centre.

The treatment is **sequential**, the pressure starts to inflate the chamber located at the end of the limb (e.g. in the legs it starts from the feet) and once the chambers are inflated, reaches the base of said limb (e.g. In the legs to the groin).

Pressoestetica® Mesis® JoySense® 2.0 - mod: JS02

Main Technical Features

Index



- **2 programs, 5 CHAMBERS**
- Adjustable treatment time
- Power: AC220V~ AC240V 50/60Hz
- 60VA consumption
- Size 340X290X300mm
- Device weight: 1,8Kg
- 24 month warranty
- EU declaration of conformity for aesthetic use

START-TIMER:
Press to start the treatment
Press after switching on to vary the time

PROGRAM SELECTION:
Slim: all full
Silhouette: lymphatic drainage massage

SWITCHING ON THE MACHINE:
Press to turn the machine on

ADJUSTING THE INTENSITY:
Turn the control knob to increase or decrease the intensity

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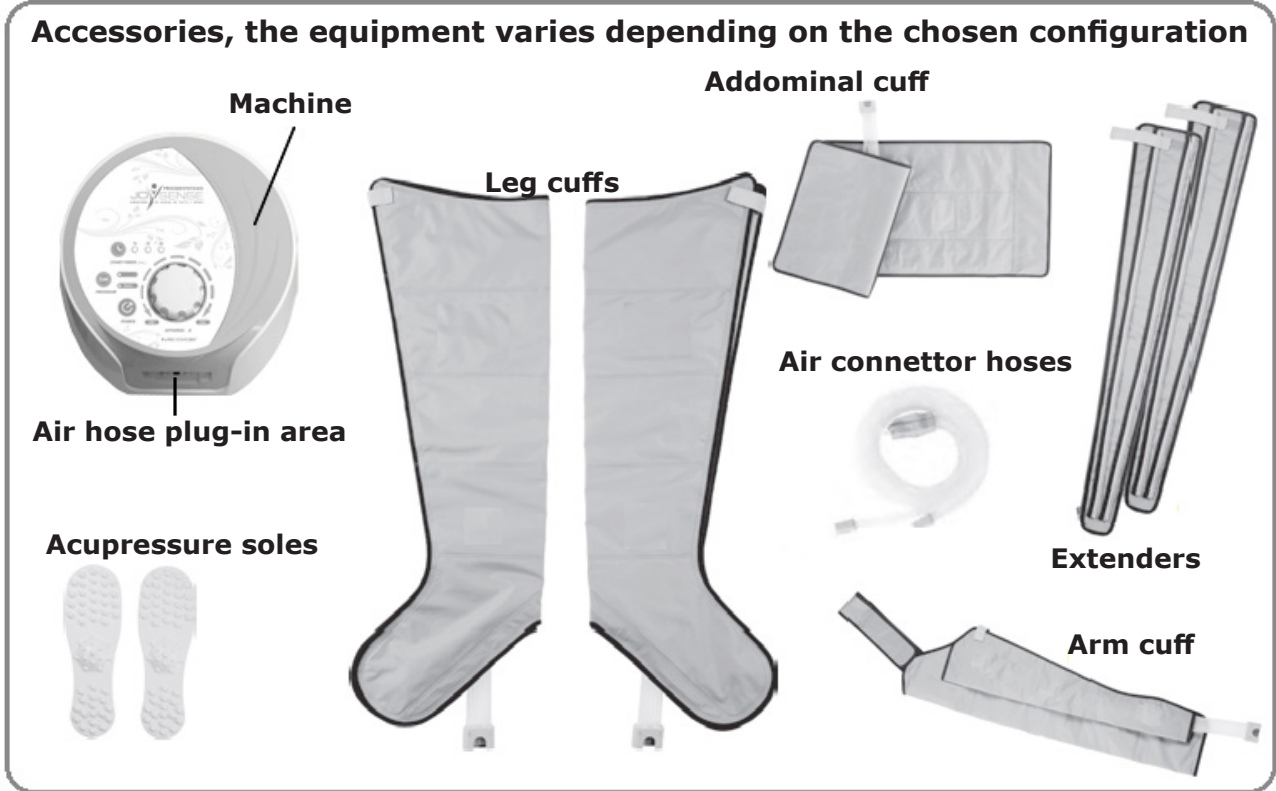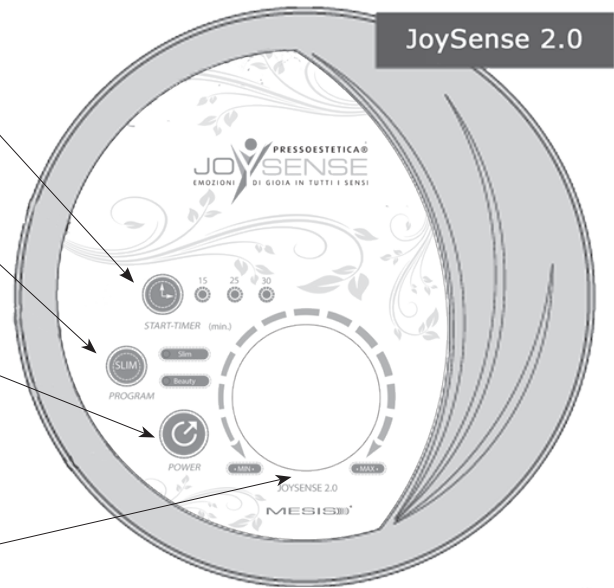- **3 programs, 5 CHAMBERS**
- Adjustable treatment time
- Power: AC220V~ AC240V 50/60Hz
- 60VA consumption
- Size 340X290X300mm
- Device weight: 1,8Kg
- 24 month warranty
- EU declaration of conformity for aesthetic use


[Index](#)

START-TIMER:
Press to start the treatment
Press after switching on to vary the time


PROGRAM SELECTION:
Slim: all full
Silhouette: lymphatic drainage massage
Beauty: deep lymphatic drainage massage

SWITCHING ON THE MACHINE:
Press to turn the machine on

ADJUSTING THE INTENSITY:
Turn the control knob to increase or decrease the intensity

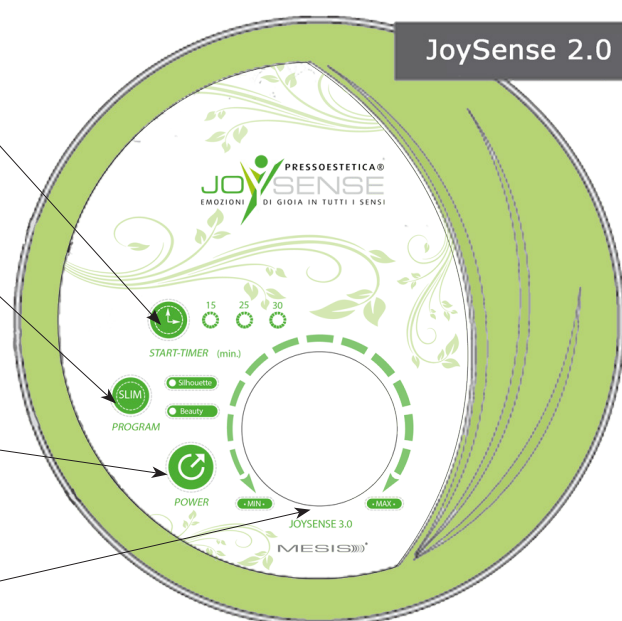
Accessories, the equipment varies depending on the chosen configuration

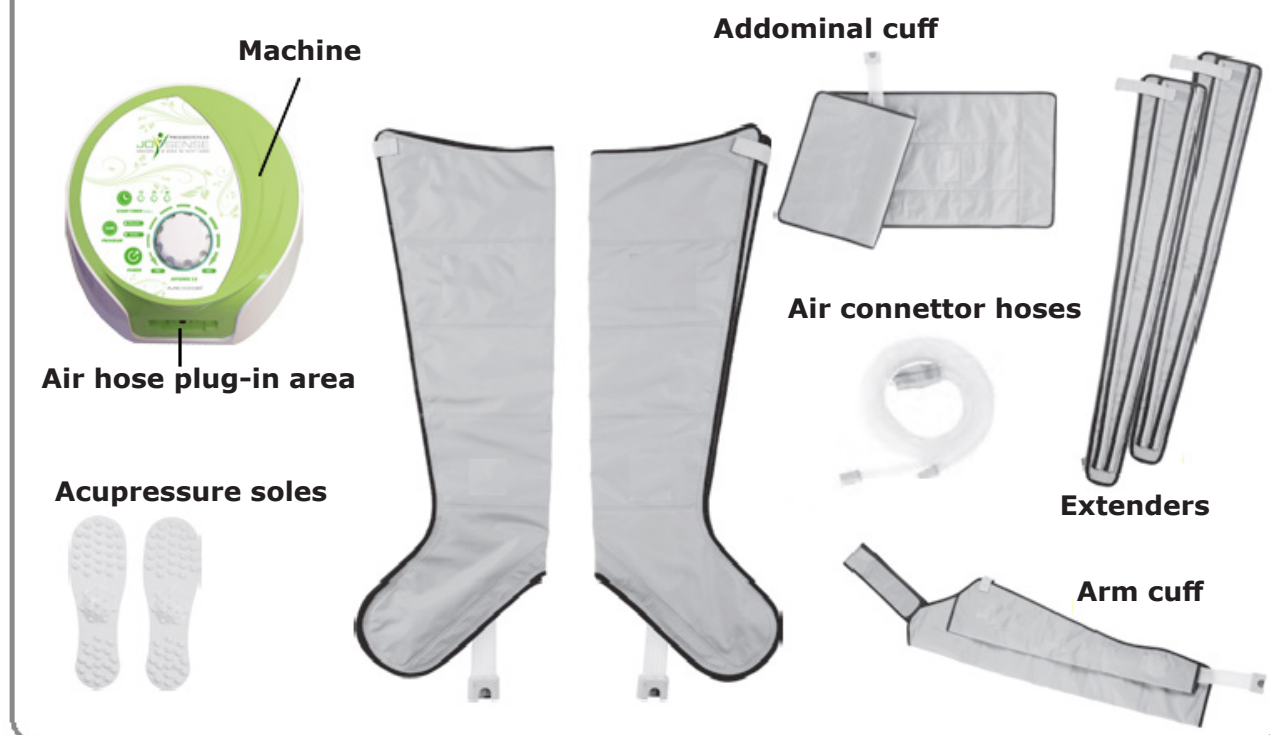
Pressoestetica® Mesis® JoySense® instructions for Use

1. Connect the power cord.

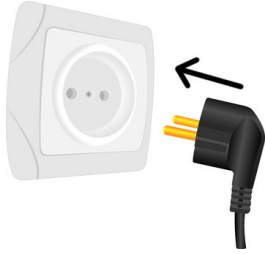
2. Connect the applicator's connector to use to the machine

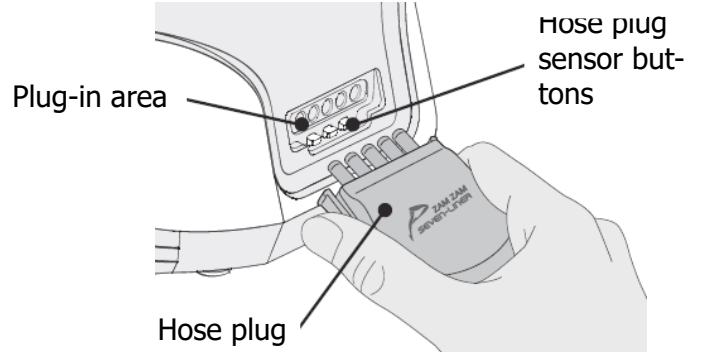3. Connect the applicators to the related connectors.

Leg cuffs



- connect the attachments of the the double connector to the leg cuffs.

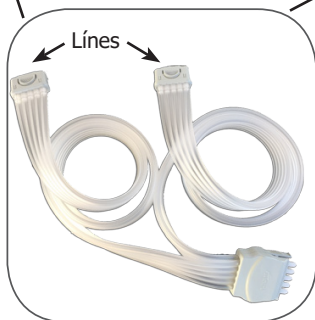Double connettor

Leg cuffs and abdominal/buttocks belt

- if you use the abdominal belt and leg cuffs, connect the double connector to the machine, after which you will have to connect one line of double connector, to the and the other one to abdominal belt. Finally, connect the two lines of the splitter to the leggings.

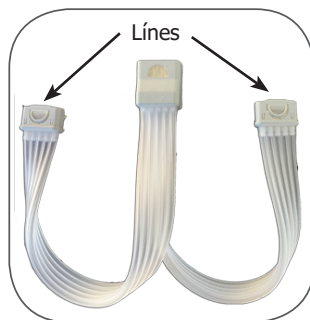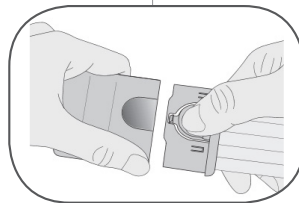Splitter

Index

Abdominal/buttocks belt

- if you use only the abdominal belt, connect the attachments of the single connector to the belt.

*Arm cuff*

- connect the attachments of the single connector to the cuff.

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- Wear the leggings and close them with the zip, place the leggings with the attachments turned inwards.
- Inflation must start at the feet and gradually reach the groin.

Abdominal/buttocks belt

- Wear the abdominal belt on the buttocks-hips-abdomen and fasten it with the Velcro.
- Inflation must start at the top part and end at the base of the abdominal belt.

Arm cuff

- Wear the armband on and fasten it with the zip.
- Inflation must start from the hand and move sequentially to the shoulder.

5. Press „POWER“ to turn the machine on.

The lamp is turned on the last used program.


6. Press „PROGRAM“ to select the program.


In the model JoySense® 3.0, when you select the program „SLIM“ remain lit the lights of both programs „SILHOUETTE“ and „BEAUTY“.

Programs JoySense® 2.0

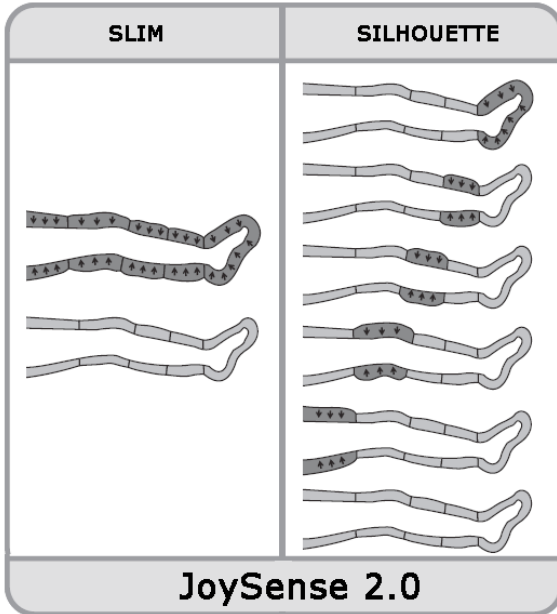Programs JoySense® 3.0

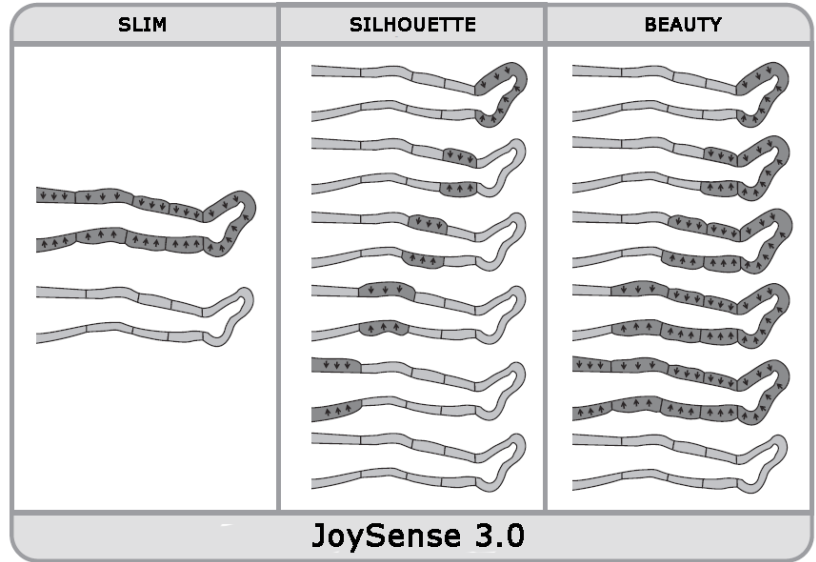
7. Press „START/TIME“ to set the duration of treatment. Press it again to change the preset time (15 - 25 - 30 mins) .

8. Turn the knob to adjust the intensity.

- To increase the intensity, turn to the right
- To decrease the intensity rotate to the left

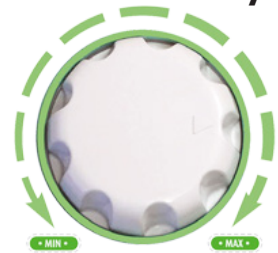N.B. When you decrease the intensity during the „SLIM“ program, inflation is stopped for about 5 seconds. In the other programs „SILHOUETTE“ and „BEAUTY“, when you decrease the intensity, the compressor automatically sets a shorter time inflation of the rooms.

Inflation starts only after have rotated the knob of the intensity adjustment.

Maintenance and Troubleshooting

The accessories cannot be repaired, as such take care of them.

Do not place the leggings, armband, abdominal belt or other accessories close to direct or indirect sources of heat.	
Do not place the equipment in damp or wet areas. Do not soak the equipment in any way, do not pour any liquid over the equipment.	
Do not wash the equipment and accessories with water or detergents. Simply wipe with a damp cloth.	
Do not place the equipment and its accessories in very cold areas or areas below zero degrees centigrade. It could cause serious damage.	
If you don't use the device for a length of time, place it back in its box.	

Environment of use and safekeeping

- **Use and maintain the device and its accessories in a clean environment and at a temperature between 10 – 40 °C, and humidity of 30 – 85%;**
- **Do not use outdoors.**

Problem found	Possible solution
The machine does not turn on	<ul style="list-style-type: none"> ● Check that it is correctly connected to the mains; ● Check that you have inserted the air tube connector to the machine; ● Check the condition of the power cord. It must be intact.
If you hear a strange hissing sound similar to an air outlet	<ul style="list-style-type: none"> ● Check that the tubes, accessories and connectors are intact; ● Check that all the connections (to the machine and the applicators) have been correctly inserted; ● Check that the air connection tubes are not bent, squashed or crushed.
The air does not flow to the applicators, there is insufficient intensity	<ul style="list-style-type: none"> ● Check that the connectors to the machine and the applicators are correctly connected; ● Check that the air connection tubes are not bent, squashed or crushed; ● Check that the chambers have inflated properly and progressively in accordance with the selected program.
The air continues to flow into the chambers and you hear an air leak in the applicators	<ul style="list-style-type: none"> ● Check for holes or damage in the applicators' air chambers. If no anomalies are detected, switch off the device, switch it back on and check that it is working properly again

 **Warranty**

- The device is guaranteed for the first user for a period of twenty four (24) months, three (3) months for the applicators and accessories subject to wear from the date of purchase against defects in materials or manufacture, provided that it is used properly and maintained in good working condition.
- The warranty is valid, if the user uses the device for professional purposes, for twelve (12) months for the machine and three (3) months for the applicators and accessories subject to wear, from the date of purchase against defects in materials or manufacture, provided that it is used properly and maintained in good working condition.

The warranty's validity is limited in the following cases:

- No warranty for accessories and materials to be deemed as "consumables"
- This Warranty is valid and effective in the country in which the Product was purchased. In the event that the product was purchased in an EU country, the warranty is still valid in all EU countries.
- Problems due to dropping the device or to improper use of the voltage.
- If the device has been opened or modified by unauthorised personnel.
- If the device has been used in unsuitable environments.
- If the device has not been used in accordance with the instructions contained in this manual.
- Warranty on the device: 2 years.
- Warranty on accessories and applicators: 3 months.

To obtain the warranty service, the user must observe the following warranty clauses:

- The products must be delivered for repair under the responsibility and at the expense of the Customer in the original packaging.
- The product warranty will be subject to the production of a tax document (sales receipt, payment receipt, sales invoice or mail-order sales document via an online shop or catalogue), attesting the product's date of purchase.
- The repair will not affect the original expiry date of the warranty and will not lead to the renewal or extensions of the warranty.
- If, during repair work, no defects are found in the product, you will still be charged the expenses relating to the technical inspection.
- The warranty is void if the defect was a result of: impact, falls, misuse or improper use of the product, use of a non-original power supply or battery charger, accidents, alterations to the product, replacement/detachment of the warranty seals and/or tampering with the product.
- The warranty does not cover damage caused during transport if unsuitable packaging is used.
- The warranty does not cover the impossibility to use the product, other incidental or consequential costs or other expenses incurred by the purchaser.

PRODUCT SUPPORT

In the event of failure or malfunction of the device, contact your dealer.

EU declaration of conformity for aesthetic use

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