



# GIMA

PROFESSIONAL MEDICAL PRODUCTS

**DEAMBULATORE IN PIEDI COMPLETAMENTE REGOLABILE - alluminio - rosso**  
**FULLY ADJUSTABLE STANDUP ROLLATOR - aluminium - red**  
**DÉAMBULATEUR DEBOUT COMPLÈTEMENT RÉGLABLE -**  
**aluminium - rouge**  
**ANDADOR DE PIE TOTALMENTE AJUSTABLE - aluminio - rojo**

**Manuale di istruzioni - Instruction Manual - Mode d'emploi - Manual de instrucciones**

È necessario segnalare qualsiasi incidente grave verificatosi in relazione al dispositivo medico da noi fornito al fabbricante e all'autorità competente dello Stato membro in cui si ha sede.

All serious accidents concerning the medical device supplied by us must be reported to the manufacturer and competent authority of the member state where your registered office is located.

Il est nécessaire de signaler tout accident grave survenu et lié au dispositif médical que nous avons livré au fabricant et à l'autorité compétente de l'état membre où on a le siège social.

Es necesario informar al fabricante y a la autoridad competente del Estado miembro en el que se encuentra la sede sobre cualquier incidente grave que haya ocurrido en relación con el producto sanitario que le hemos suministrado.

**Gima 43172**



**Foshan HCT Medical Equipment Co., Ltd**

No.501 of 19 Building, Daxing Liandong U Valley, No. 4, Xiangshun Road, Danzao Town, Nanhai District, Foshan, Guangdong, China

Made in China

**REF**

**HCT-9291D**

**EC REP**

SUNGO Europe B.V.

Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel, The Netherlands

**MD**



**Gima S.p.A**

Via Marconi, 1 - 20060 Gessate (MI) - Italy  
gima@gimaitaly.com - export@gimaitaly.com  
www.gimaitaly.com



## WARNINGS

Before carrying out any operation, visually check the integrity of the device. Where the product presents anomalies, immediately contact the dealer from whom it was purchased. We assume no responsibility for the misuse of the device.

## DESCRIPTION UPRIGHT WALKER

HCT-9291D upright walker features compact folding design and aluminum stable frame and fully adjustable armrest, aimed at helping users to stand upright and Look ahead with better support.

*The Upright Walker requires assembly before use. Please follow the instructions below:*

## MATERIAL COMPONENTS

<b>Frame</b>	Aluminum
<b>Wheels</b>	PVC
<b>Handrail</b>	Foam
<b>Handrest</b>	Padded PU material
<b>Seat</b>	Nylon and polyester

## PREPARING FOR USE

The choice of the type of aid and the necessary settings for use must be carried out by medical and technical staff qualified, depending on the pathology and user terms and conditions.

1. To ensure your safety in using the the upright walker, the following safety information and all instructions must be followed.
2. Check whether all wheels are tight before use.
3. Check whether the height adjustment button is firm.
4. Make sure the walker is fully deployed. Pay attention to the surrounding area when using a walker.
5. The backrest is intended to provide back support while seated. The backrest is not intended to support full body weight. Avoid Leaning while sitting as this may result in fall.
6. The brakes MUST be in the locked position BEFORE sitting on the seat.
7. Brakes are to be used to keep the wheels from rolling. However, it may not prevent the walker from moving if it is pushed.
8. DO NOT self-propel the walker while seated. Walker are intended to use as a walking aid for individual only and are NOT TO BE USED AS A WHEELCHAIR.
9. DO NOT use the walker to walk backwards, down gradients or to climb stairs, curbs or to go over big obstacles. Serious risk of fall or injury may occur.
10. DO NOT use the seat to transport people or Large objects.
11. DO NOT sit on the seat when the walker is on an incline.
12. DO NOT push the walker while brakes are enabled.
13. DO NOT put your hand into the folding place of the walker to avoid being pinched.
14. DO NOT Let yourself fall on the seat of the walker, but always sit slowly on the seat in a controlled manner.

## COMPONENT



1. **Handgrip**
2. **Backrest**
3. **Break wire**
4. **Shopping nylon bag**
5. **10 inch front wheels**
6. **Quick-release screw**
7. **8 inch rear wheels**
8. **Obstacle Aid Pedal**
9. **Brake cover**
10. **Cane holder**
11. **Joint part for folding design**
12. **Padded Seat**
13. **Push button for height adjustment**
14. **Foam handrail**
15. **Knob to adjust backrest height**
16. **Knob to adjust handle Length**
17. **Arm holder**

## TECHNICAL SPECIFICATIONS

Width between the handles	19.7 in
Overall width	24.6 in
Overall depth	29.5 in
Seat area	9.8*18 in
Seat height	21.6 in
Adjustable handle height	40.3-45.5in
Max. user weight	300Lb (136kg)
Weight	22.4 lb
Front wheels	10in
Rear wheels	8 in

## SAFETY INSTRUCTIONS

The Upright Walker requires assembly before use. Please follow the instructions below:

1. Identify the partially assembled base piece, unfold it, and ensure the seat is flat and Locked into place.



2. Attach the 2 front wheels to the seat inserting and Lock it into position.



3. Locate the 2 handles which comes with a backrest and insert them into the two tubes at two sides of the walker. Make sure it is Locked into the position.




4. Insert the 2 hand grips into the designated slot within the arm bars. Choose the desired Length by Lining up the scale number where you want the Length to be and secure with the knob by screwing it tight.



5. Put handrail into the tube and Lock into the position.



6. Pull out two roll tubes at the front of the seat, then hang the bag through the two straps of the storage bag.

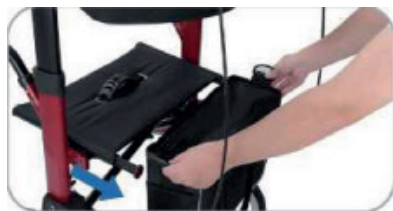

## HOW TO ADJUST ARMHOLDER HEIGHT

There is a push button at each tube where the handle height can be adjusted up or down by pressing the button and snapping into place when desired height is reached.



## OPERATING THE HAND BRAKES

Pull in the hand brake to slow or stop the walker. Push out the hand brakes until they "click" to lock the walker in position.

Please attention: Brakes must be applied in the lock position when seated. Pre-check the Lock by pushing forward on the upright.



## FOLDING THE UPRIGHT WALKER

To fold the upright walker, simply pull up the strap at the middle of the seat. Then pull back forward down the upright walker into flat.

There are two screws at each sides of handle bar, backrest height can be adjusted up or down by Loosen the screws and tighten them when desired height is reached.



## HOW TO ADJUST THE HEIGHT OF THE BACKREST



## BRAKE ADJUSTMENT:

Turn the adjusting screw anti-clockwise to increase the tension then Lock the adjustment Locking nut clockwise whilst holding the adjuster screw stationary.

## MAINTENANCE

Over a period of time the brake cables may stretch, and if no further adjustments can be made the cables should be replaced. Seek the help of your dealer when replacing brake cables. The castors should be checked that they are secure each time before using the walker. Also, check the brakes operate correctly and the handles and hand knobs are secure. Should any part of the walker appear Loose or insecure discontinue use until it is inspected by an authorized dealer.

The following shows the maintenance schedule:

### WEEKLY

- Check brake for malfunction
- Check all wheels and castors for Loose bolts and nuts

### PERIODICALLY

- Check handgrips for wear
- Check frame bolts are secure
- Check front and rear axle bolts are secure
- Check tyres for wear or cracks

### YEARLY

Once a year take your walker to an authorised repairer for a full check, clean, service and Lubrication of all moving parts

Wheels Hand grips Brake Cables Brake Lever Bearings Knobs

## CLEANING

Your walker can be cleaned with a damp cloth or sponge and warm soapy water. Domestic anti-bacterial surface cleaners may be used if required. Care should be taken not to leave any cleaning residues on the handgrips or seat which may Leave them slippery when used. Do not use scouring pads or scrapers as this may damage surface coatings.

















### WARNINGS!

1. Maximum working Limit 300Lb (136kg).
2. Do NOT stand on this product.
3. Ensure both push handles are adjusted to the same height.
4. Seek advice from a trained health care professional or approved Freedom Healthcare dealer for height adjust merit of this product.
5. Always maintain a safe distance between your walker and other pedestrians or objects.
6. Some components may cause harm or damage when they come into contact with other objects.
7. Do not attempt to replace the wheels or castors of your walker. If they are faulty or require a replacement please contact an authorized repairer.

8. Ensure that your walker is opened fully (folding scissor arm is straight) before using it.
9. Ensure that your fingers are clear from the scissor arm when folding.

### UPRIGHT WALKER USER MANUAL

We warrant the included accessories against defects in materials and workmanship for one year from the date of original retail purchase. We do not warrant against damage caused by accident or abuse. If you have any questions, please contact us.

	Caution: read instructions (warnings) carefully		Medical Device compliant with Regulation (EU) 2017/745		Product code
	Manufacturer		Non-sterile		Consult instructions for use
	Keep away from sunlight		Keep in a cool, dry place		Medical Device
	Importer		Authorized representative in the European community		Date of manufacture
	Lot number		Unique Device identifier		

### GIMA WARRANTY TERMS

The Gima 12-month standard B2B warranty applies.