457

Pulse Oximetry

The world's #1 oximeter

9550 Myx²1

All-in-one fingertip pulse oximeter

Vital Signs Within Reach.™

Product Highlights

- ✓ Simple Automatically turns on for fast, easy spot-checks
- ✓ Superior Accuracy Trusted by clinicians worldwide
- ✓ Efficient New turn off mechanism allows for 2500 spot-checks or 21 hours of continuous use
- ✓ Cost Effective No additional sensors needed
- ✓ Enhanced Design Enclosed spring; smaller; wider viewing angle; easy battery access

Nonin, the world leader in fingertip oximetry, has designed years of clinical experience into the Onyx II. The electronics and sensor are incorporated into the world's smallest fingertip oximeter, providing a cost effective solution for spot-checks and short term monitoring. Onyx II – not just a new look. Superior performance and unparalleled reliability put vital signs within reach anywhere, anytime.



Vital Signs Within Reach

The Onyx II is ideal for use in any situation where a fast accurate reading of blood oxygen saturation and pulse rate is needed. Never search for a pulse oximeter or sensor again.





Proven Performance/ Simple Operation




✓ Small, Lightweight - Carry the Onyx II anywhere on a lanyard, in a carrying case or in your pocket. It weighs less than 2 oz.

✓ Large LED display - The widest viewing angle available for optimum visibility of patient readings, day or night.

✓ Battery Door - Enhanced design for convenient battery access.

✓ Pulse Quality Indicator - Brilliant, easy-to-see LED provides a quick and

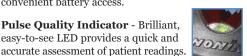





✓ Energy Efficient - A new, turn-off mechanism allows two AAA batteries to provide approximately 2500 spot checks or 21 hours of continuous operation.





Pediatric to Adult Finger Sizes - Accommodates a wide range of finger thicknesses, 0.3-1.0 inch (8 mm to 26 mm).

✓ Durable - So rugged and reliable, it can be used in the most demanding environments.

Everything You Need





Standard Package Items

Please contact your authorized distributor for more information.

Specifications

Oxygen Saturation Range (%SpO₂) 0% to 100%

Pulse Rate Range 18 to 321 pulses per minute

Displays

Numeric Displays Two 3-digit LEDs, 7 segment LED (light emitting diodes) Pulse Strength Tricolor LED segments

Red 660 nanometers Infrared 910 nanometers

Measurement Wavelengths

Blood Oxygen Saturation (%SpO $_2$ ±1 SD)^a 70 - 100% ± 2 digits Pulse Rate 20 - 250 bpm ±3%

Temperature

Operating +32° to +122°F (+0° to +40°C) Storage/Transportation -22° to +122°F (-30° to +50°C)

Humidity

Operating 10% to 90% noncondensing Storage/Transportation 10% to 95% noncondensing

Operating Altitude Up to 40,000 feet Hyperbaric Pressure Up to 4 atmospheres

Power Requirements Two 1.5 volt AAA-size batteries

Battery Life

Approximately 21 hours of continuous operation Operating Storage

4 years, minimum

Dimensions 1.31 x 2.22 x 1.27 in. (3.32 x 5.63 x 3.22 cm)

Weight 1.89 oz (53.5 g) with batteries

Warranty 2 years

Classifications per IEC 60601-1

Degree of Protection Type BF-Applied Part

S.D. (Standard Deviation) is a statistical measure; up to 32% of the readings may fall

^bUsing two AAA-size alkaline batteries; based on approximately 2500 spot checks calculated at 30 seconds per spot check. 'Using two AAA-size alkaline batteries.

Quality Systems are registerd to ISO 13485:2003 and EC certificates per Annex II, clause 3 of EC Directive No. 93/42/EEC concerning medical devices.



Specifications are subject to change without notice.

Your Partner of Choice

With over 18 years of experience and dedication to the design and support of noninvasive monitoring devices, Nonin has helped many medical professionals meet clinical and economic objectives. Trusted by clinicians worldwide, countless tests are performed every day on thousands of Nonin pulse oximeters in more than 125

countries. Nonin's dedication to technological leadership, precision manufacturing and uncompromised customer support ensures quality products and service you expect.



Nonin Medical, Inc. 13700 1st Avenue North Plymouth, MN 55441-5443, USA Phone: +1.763.553.9968 Toll Free: 1.800.356.8874

Fax: +1.763.553.7807 info@nonin.com